


# Domus 4

## Exceptional comfort with advanced digital control for pressure redistribution system

Experience unparalleled comfort with Domus 4, our advanced air mattress equipped with a digital control panel. This cutting-edge design provides precise control over air levels, delivering consistent support to prevent bottoming out, and ensure coccyx pressure relief even while seated. Additionally, the function offers a stable surface for caregivers to perform tasks easily, enhancing overall functionality and convenience.


Domus 4 ensures a tailored sleep experience, enhancing immersion and envelopment for a soothing sleep environment. Plus, the mattress operates quietly, minimizing disturbance and contributing to overall satisfaction.<sup>1</sup>

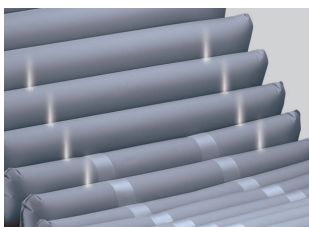


**Microclimate Management**  
Micro low air loss effectively wicks away moisture from the surface, leading to a more comfortable experience.



**Protect Vulnerable Heels**  
The individual air cell can be easily removed to achieve zero pressure, preventing heel pressure injuries.

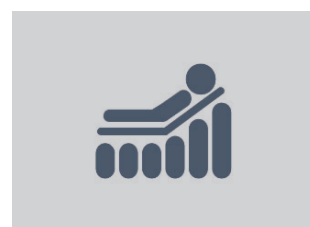

**Prevent Bottoming Out while Sitting**  
**Seat Inflate Mode**  
manually manages the inflation of the air cells to improve stability and comfort during sitting.

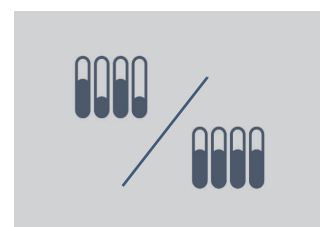


**Therapy Modes to Meet Diverse Needs**  
Provide Static mode for enhanced patient immersion, Alternate mode for allowing pressure redistribution, and Max Firm mode for nursing procedures.

# Domus 4

Personalized comfort and pressure injury prevention at your fingertips



**CPR Knob**  
Quickly deflate the mattress for emergencies.

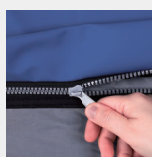

**Decreases the Risks of Infection.**  
The top cover is made with welded seams (high frequency induction hardening ®) to avoid fluid ingestion.




**Special Design of Air Cells**  
Cell-on-cell design for continuous pressure support, even during power failures.



**Storage Strap**  
Side buckle under the base cover can be used to store rolled mattress very easily.



**Quick Connector with Cap**  
Mattress can be transformed from a powered support surface to a sealed static air mattress for transportation need.

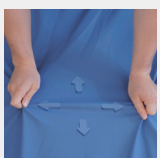

**4 Way Stretch Top Cover**  
Offering low shear, low friction, vapour permeability, and moisture protection, coupled with optimal immersion and envelopment.



Water Resistance




Flame Resistance




Biocompatibility




Vapour Permeable

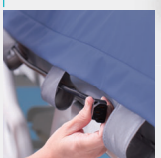

**Cable Management**  
Five loops on the side of mattress for managing cable to avoid tripping hazard.

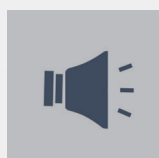

**Prioritize Patient Safety**  
Visual and audible signals are used to alert caregivers in case of any failure.

Specifications		Domus 4				
<b>Pump</b>		Dimension	29 x 18.5 x 12.6 cm (11.4 x 7.3 x 5 in)			
		Weight	2.3 kg ( 5.1 lb )			
		Case Material	Fire Retardant ABS			
		Supply voltage	220 - 240 V / 50 - 60 Hz			
		Therapy mode	Max firm / Alternate / Static / Seat inflate			
		Noise Level	< 30dB			
<b>Mattress</b>		Dimension	Cells	Length	Width	Height
			17-20	200cm ( 78.7 in )	85-90cm ( 33.5-35.4 in )	12.7-20.3cm ( 5-8 in )
		Type	5" + 3" foam or 8" Replacement			
		Weight	5.6 - 6.9 kg ( 12.3 - 15.2 lb )			
		Top cover material	Nylon / PU, Poly / PU			
		Cell material	Nylon / TPU			
		Maximum patient weight	180 - 250 kg ( 396 - 550 lb )			

Pump: water resistant standards (IP21); Mattress: flame retardant standards (EN597-1, EN597-2), RoHS, WEEE

